

# Desara<sup>®</sup>TV

TRANSVAGINAL SLING

**2.7mm**  
Flexes 1 cm less than the 2.7mm competitor product<sup>1</sup>.

**3.0mm**  
Flexes 2.4 cm less than the 2.7 mm competitor product<sup>1</sup>.


Two Trocar  
Diameters

Tighter trocar curve closely contours the pubic bone while maintaining distance from critical anatomy

85° Curve

Enclosed Dilator Tips

Smoother tissue passage and reduced tissue drag, resulting in the lowest insertion force when compared to competitors<sup>1</sup>.

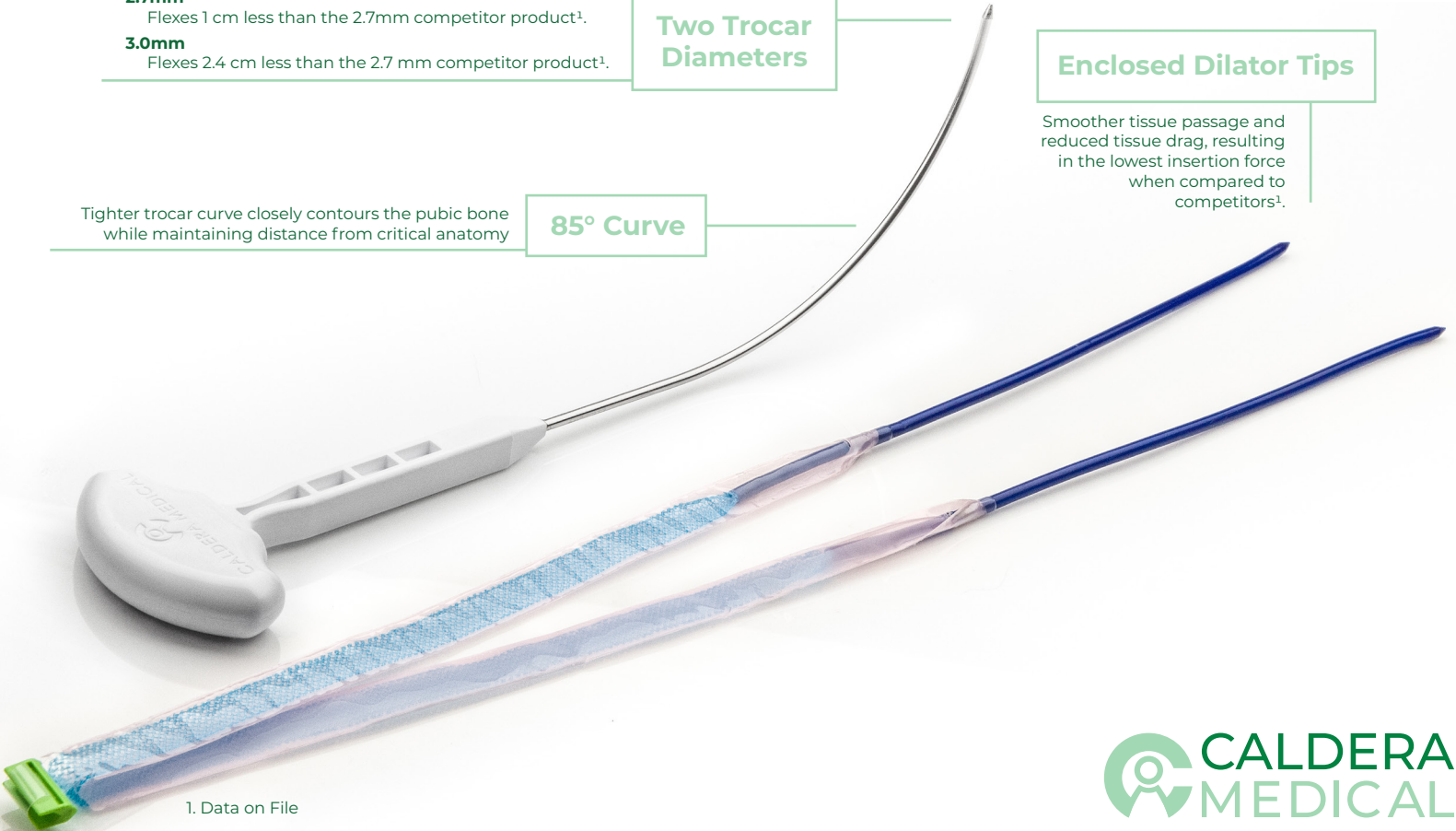

1. Data on File

## Disposable Instrument

5mm increase in trocar height for higher BMI women

## Simplified Dilator Tube Removal

Extension of dilator tube facilitates easy trocar removal



## EZ Tab

1. CUT midline sleeve at centering tab
2. REMOVE the highly visible one-piece EZ tab
3. PULL dilator tubes to remove sleeve

# Desara®TVez

TRANSVAGINAL SLING

Desara® mesh is engineered to lay flat when placed under tension and does not fray or curl<sup>1</sup>. Offering stable and controlled support.

Blue mesh allows for better physician visualization, compared to clear slings.

Type I large pore mono-filament polypropylene designed with approximate pore size of 1000 microns to promote tissue in-growth and allow for macrophage movement.

## ORDER INFORMATION

ORDER NUMBER	PRODUCT DESCRIPTION	QUANTITY
CAL-DSO1BTV27	DESARA TVez Blue System with 2.7mm Introducer	Order in Multiples of 5
CAL-DSO1BTV30	DESARA TVez Blue System with 3.0mm Introducer	Order in Multiples of 5

10-462 Rev A  
1. Data on File

